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User manual
Digital Refractometer
Ref. 4817

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1. Instrument structure

2. Instrument measurement


- Before the testing
- Power on
- Power off
- Operating tests
- Zero point calibration

3. Note

4. Technical parameters

5. Common trouble shooting

1. Instrument structure

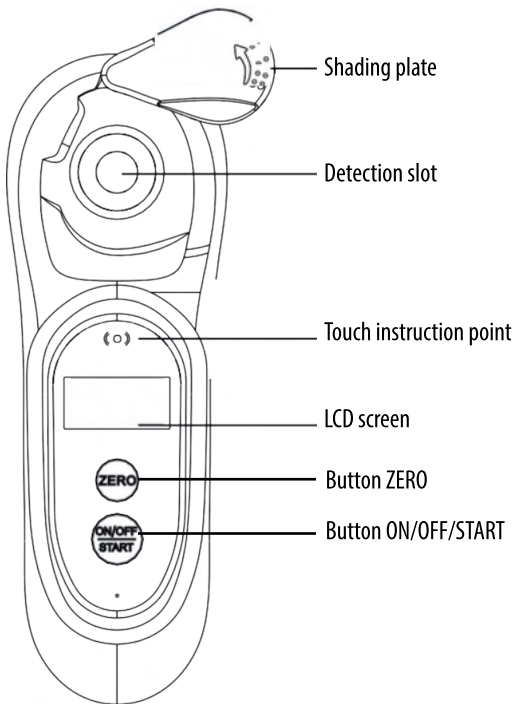2. Instrument measurement

Before the testing :

The instrument and liquid test sample should be placed in same environment temperature. (The temperature difference between the environment, the instrument, and the liquid sample must not exceed 2°C). Otherwise it will affect the accuracy of the testing.

Shake the liquid sample thoroughly before testing..

If the instrument has not been used for a long time, resetting the "Zero point" before using.

Check the battery level. Low battery power may affect measurement accuracy.

Power on :

Press and hold the "ON/OFF/START" for one second until the screen display.

Power off :

Press and hold the "ON/OFF/START" for at least one second until the screen display "Power off".

Operating steps :

Open the shading plate counterclockwise.

Use pure water and a soft cloth to clean the detection area, then remove the residual water.

Add 3~5 drops of liquid test sample to the detection slot.

Close the shading plate clockwise.

Press the "ON/OFF/START" shortly, or touch the "touch instruction point".

The screen shows the value of BRIX and Value of Refractive index.

Zero point calibration:

Generally the instrument can be used directly. Zero point calibration need to be performed when the instrument is long time unused or there's a zero offset.

The method is as follows:

Use the pure water to clean the detection slot with soft cloth repeatedly.

Add 3~5 drops of pure water to the detection slot.

Press and hold the "ZERO" until the screen shows "zeroing", release the button and start.

When the Screen shows "Zero Success", it means the Zero point calibration is finished. When the screen shows "Zero Failure", it means the Zero point calibration is failed.

Press the "ON/OFF/START". If the displayed value is 0.0%, the zero-point calibration is successful.

3. Note

Make sure the clean pure water is used for Zero Point Calibration.

Make sure to clean the detection slot after testing.

The instrument should be placed in the environment of dry and normal temperature. Avoid exposing the instrument to high temperatures or direct sunlight to protect measurement accuracy and the OLED display.

Don't wash the instrument immersion into the water to avoid the leaking.

Digital Refractometer is a precision optical instrument, please using and maintaining carefully and softly. Avoid the damage to the prism in detection slot from cut and bump.

Avoid sample evaporation, as it may affect the measurement results, please operate fast.

Low battery power may cause incorrect measurement data.

4. Technical parameters

Parameter	Value
Measuring range	0.0–95 Brix (%)
Resolution	0.1 Brix (%) / 0.1 °C
Accuracy	±0.3 Brix (%) / 1 °C
Temperature compensation	10–80 °C
Operating temperature	10–40 °C
Sample volume	~0.5 mL (3–5 drops)
Measure time	~3 s
Power supply	2 x AAA
Dimension	154 × 52 × 44 mm
Weight	140 g

5. Common trouble shooting

- Concentration is too high or too low :

Situation	Paraphrase	Cause	To do
Concentration is too high or too low.	Outside the measurement range	The sugar content of the tested sample is higher than the detection range	Dilute the testing sample.
		Influence of external light	Close the shading plate properly
		The tested liquid contains particles	Filter out the particles

Situation	Cause	To do
Concentration is too high or too low.	Liquid temperature is lower or higher than the instrument range	The temperature of the liquid should be the same as the temperature of the instrument and the environment
Pas de mesure	Missing sample solution	Add the test solution to be tested.

- Zero failed :

Cause	To do
Detection slot is not clean	Cleaning the detection slot.
The calibration solution is not correct	Change the calibration solution.

- Test data is incorrect :

Cause	To do
The surface of the prism is uncleaned	Wipe or clean prism
Improper operation	Please follow the instructions
Photometric effect	Operation after covering the shading plate
Operation at high temperature	Operation at normal temperature.
Low battery power	The power is displayed as empty, which will affect the test data.

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Cause	To do
Misjudgment of product quality due to improper operation	Please follow the instructions